


Application Note

Aromatic Compounds

Baseline - MOCON, Inc.™

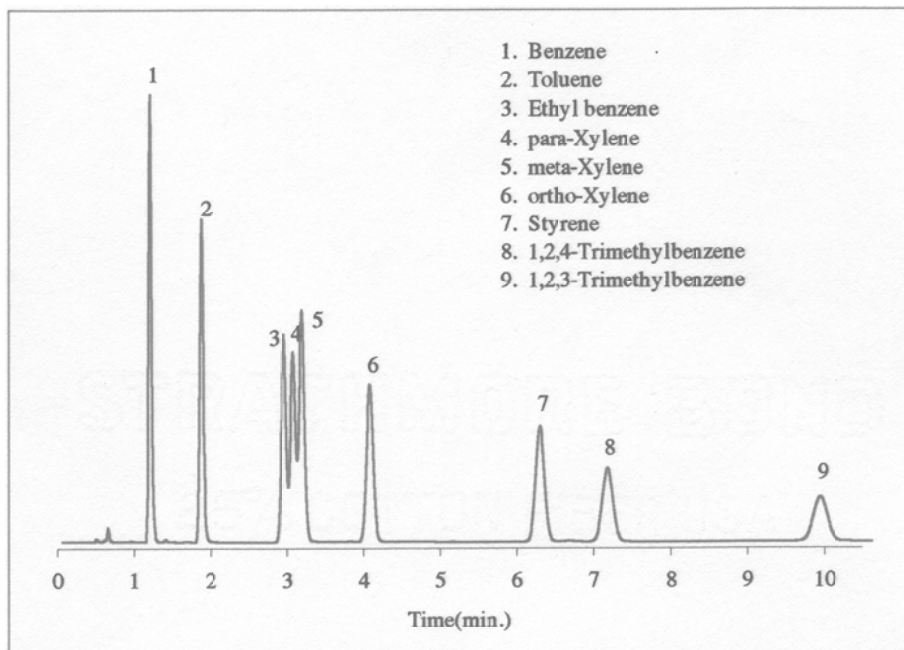

Chromatogram

Analysis

Detector: FID

Represented by:

Column: DB Wax, FSOT, .53mm ID x 30m


Oven Temperature: 60°C

Carrier: Helium, 25cc/min typical

Fuel: Hydrogen, 40cc/min typical

Combustion Air: 200cc/min typical

*Baseline, the reference point
from which all things are measured.*

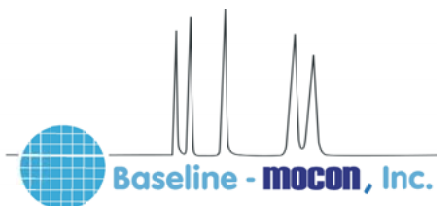

19661 Highway 36, P.O. Box 649, Lyons, CO 80540
In the continental United States, phone 800.321.4665, or fax 800.848.6464, toll free.
Worldwide, phone 303.823.6661 or fax 303.823.5151
• URL: www.baselineindustries.com • E-mail: sales@baselineindustries.com